



Apply and release brake

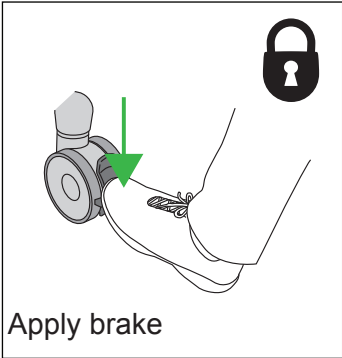


Apply brake

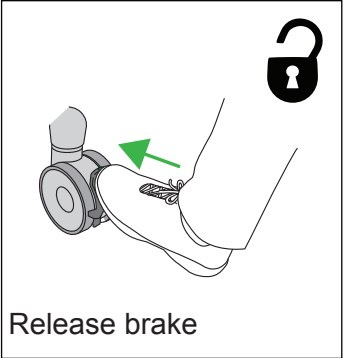


Release brake

Straight steering device

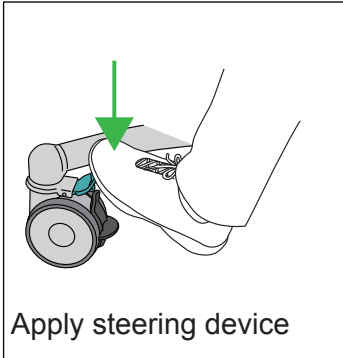


Apply steering device

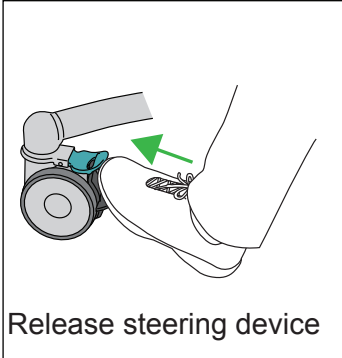


Release steering device


Control panel

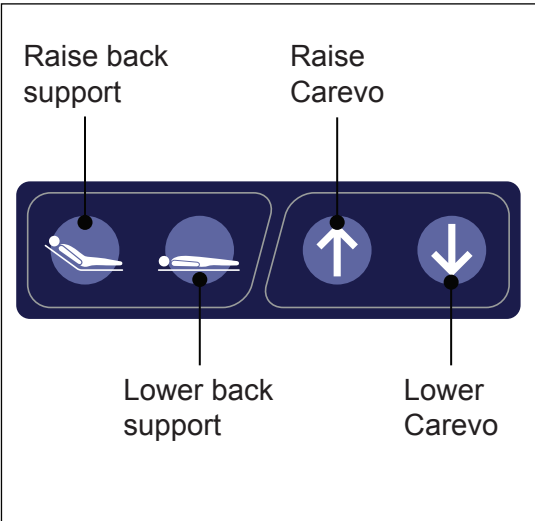

Side support

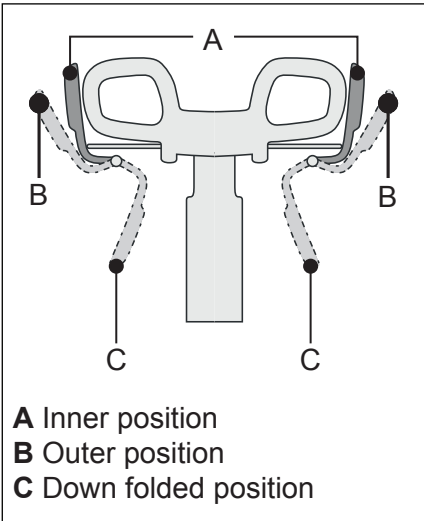


Emergency stop button

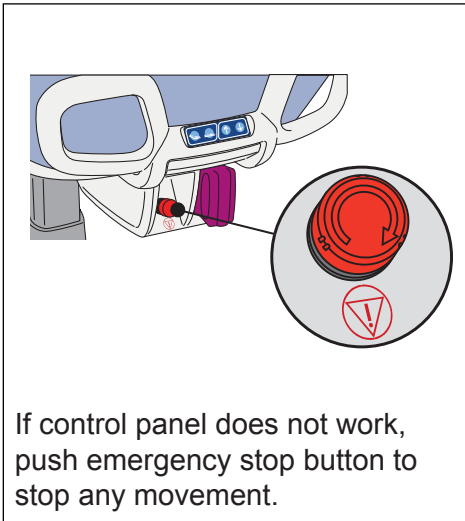


If control panel does not work, push emergency stop button to stop any movement.

Battery insert

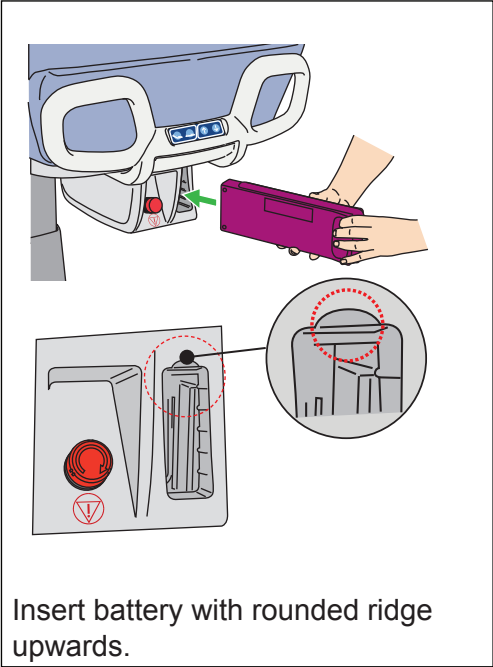


Insert battery with rounded ridge upwards.

Battery remove

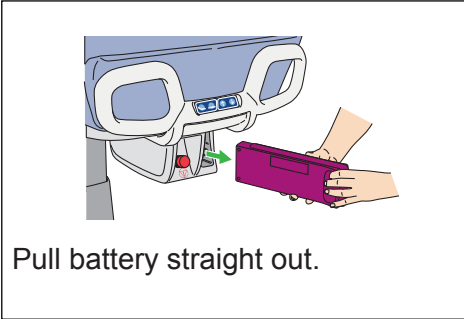


Pull battery straight out.

Comfort handles

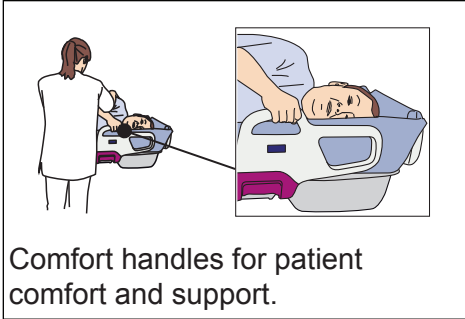


Comfort handles for patient comfort and support.

Flexi zone

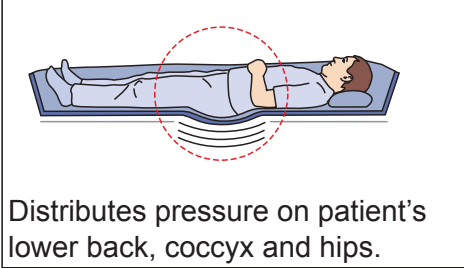

Distributes pressure on patient's lower back, coccyx and hips.

Flexi zone position

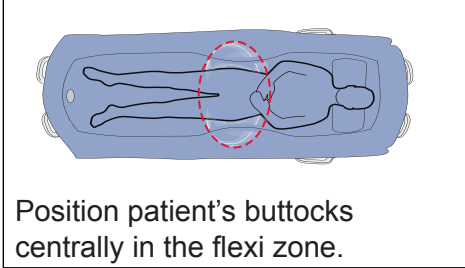

Position patient's buttocks centrally in the flexi zone.