



FD1800, FD1805

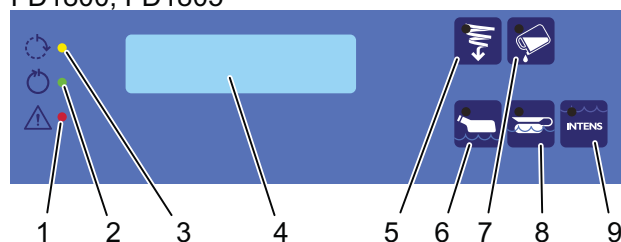

FD1810

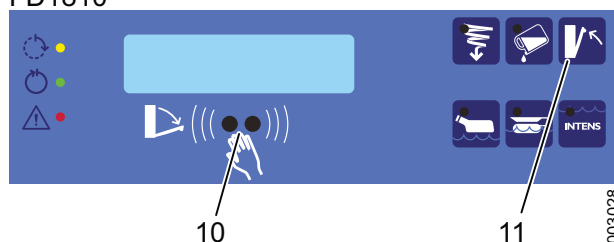

	Function
1	Red: Error indication. The display shows an error message.
2	Green: The program is complete. The light goes out when the door is opened and for a short time the display shows a program summary in the window, for example maximum temperature.
3	Flashing yellow: The program can be terminated. Yellow: A program is running. The display shows the current status, for example the cleaning phase and temperature.
4	Display. A service technician can enter the information to be shown in the display.
5	Button for selecting rim flushing.
6	Button for choosing economy program for lightly soiled goods.
7	Button for selecting cleaning agent (option).
8	Button for choosing normal program for normally soiled goods.
9	Button for choosing intensive program for heavily soiled goods.
10	IR sensor to open/close the door and to start the preselected program (option).
11	Button for opening/closing door (option).

Code	Description	Measure
U1	Lime descaler low level	Top up descaler agent.
U2	Process agent low level	Top up process agent (option).
U3	Rinse-aid low level	Top up rinse-aid (option).
U6	Inspection request	Contact authorized service personnel.
U7	Warning of high goods temperature	Let the goods cool down.
U8	External communication	External connection interrupted.
U9	Inactivity warning	The machine has been inactive for 72 hours. A program with an empty chamber must be run before using the machine to clean goods.
H1	Descaler empty	Top up descaler agent.
H3	Rinsing agent empty	Top up rinsing agent.
H6	Start attempt with door open	Close the door (the machine makes a new start attempt). If the error recurs, contact authorized service personnel.
H7	Battery low	Contact authorized service personnel.

For type F error codes, see Instructions for Use.