


INSTRUCTIONS FOR USE

# Wall Mounted Battery charger

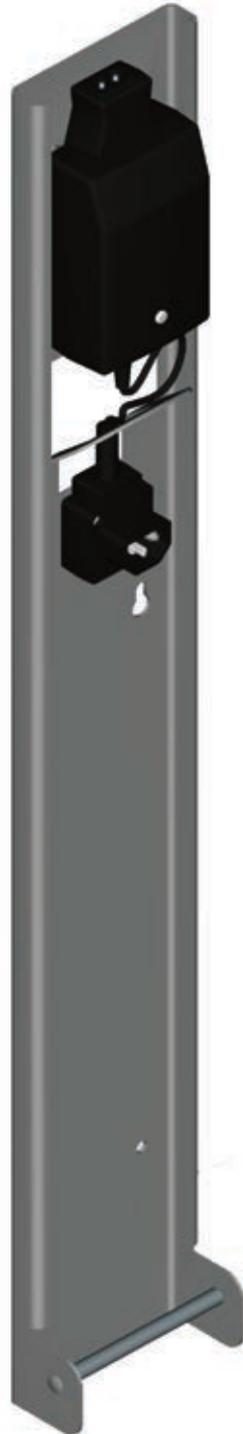

# Table of Contents

---

<b>General Information</b> .....	<b>3</b>
Foreword .....	3
Service and Support .....	3
Manufacturer Information .....	3
Definitions Used in this Manual .....	3
Intended Use .....	3
Operational Life .....	3
Product Identification .....	3
<b>Safety Instructions</b> .....	<b>4</b>
Battery and Battery Charger Safety Practices .....	4
<b>Installation Instructions</b> .....	<b>5</b>
Preparation .....	5
Installation .....	5
<b>How to Use the Charger</b> .....	<b>5</b>
Charging your Battery .....	5
<b>Troubleshooting</b> .....	<b>6</b>
<b>Technical Specifications</b> .....	<b>6</b>
Technical Data .....	6
Dimensions .....	6
<b>Electromagnetic Compatibility</b> .....	<b>8</b>
Electromagnetic Compliance.....	8
Electromagnetic Immunity.....	9

Design Policy and Copyright

® and ™ are trademarks belonging to the Arjo group of companies.

© Arjo 2019.

As our policy is one of continuous improvement, we reserve the right to modify designs without prior notice.

The content of this publication may not be copied either whole or in part without the consent of Arjo.

# General Information

Thank you for purchasing the Wall Mounted Battery Charger.

Your charger is part of a series of quality products made for hospitals, nursing homes and other health care uses.

Our goal is to serve your needs and provide you with the best products available. We also offer trainings that will bring your staff maximum benefit from every Arjo product.

Please contact us if you have any questions about usage or maintenance of your Arjo product.

## Foreword

Please read this manual in its entirety before using your charger. The information in this manual is crucial for the proper usage and maintenance of the charger. It will help to protect your product and make sure it performs to your satisfaction. This manual must be read and understood as it contains safety related topics that may help you to prevent injuries.

**WARNING: Use only this product with Arjo designated parts to avoid any lack of performance or injuries.**

If a serious incident occurs in relation to this medical device, affecting the user, or the patient then the user or patient should report the serious incident to the medical device manufacturer or the distributor. In the European Union, the user should also report the serious incident to the Competent Authority in the member state where they are located.

Modifications to any Arjo product may affect its safety. Arjo will not be held responsible for any incidents that occur as a result of any unauthorized modification to its products.

## Service and Support

For more information, contact your Arjo agent which can offer support to maximize the long-term reliability, safety and value of the product. Contact your local agent for product replacement. This product does not require maintenance, and have no parts that can be replaced.

You can get additional copies of this manual from your local agent. When ordering, include the *Instructions for Use* product number (see front page) and product part number.

## Manufacturer Information

This product is manufactured by:

ArjoHuntleigh AB  
Hans Michelsensgatan 10  
211 20 Malmö,  
SWEDEN

☎ : +46 (0) 10-335 45 00

📠 : +46 (0) 413-138 76

🌐 : www.arjo.com

## Definitions Used in this Manual

### WARNING:

Means: Failure to understand and follow these instructions may result in injury to yourself and others.

### CAUTION:

Means: Failure to follow these instructions may cause damage to the product.

### NOTE:

Means: This is important information regarding the correct use of the equipment.

## Intended Use

The charger is intended for charging Arjo batteries according to instructions herein. All other uses must be avoided.

The charger is for indoor use only, under normal conditions (0 °C to 40 °C [32 °F to 104°F], 90 % RH max).

## Operational Life

The charger has been made and tested for a life span of ten (10) years based on a usage rate of two (2) complete recharge (8 hours) per day. The life span corresponds to the safe operating time period before the product replacement. Aging of the unit and frequency of use are factors that have an impact on the product's life span.

## Product Identification

The unit's identification number (specification, model, serial number) appears on a silver label attached to the charger (see Fig. 4).

# General Information

## Symbols Used

Symbols	Key to Symbols
	This symbol points out the manufacturer's name and address. May also point out manufacturing date.
	CE marking indicating conformity with European Community harmonised legislation.
	This symbol means that the product must be kept dry.
	This symbol points out the manufacturer's catalogue number.
	This symbol points out the manufacturer's lot number.
	This symbol means that you must refer to the <i>Instructions for Use</i> .
	Waste Electrical and Electronic Equipment (WEEE) - do not dispose of this product in general household or commercial waste.
	The device is meant to be operated indoors only.

Symbols	Key to Symbols
	Class II electrical equipment. This term refers to electrical equipment in which protection against electric shock does not rely on basic insulation only.
	This symbol means the product is certified by the Slovenian Institute of Quality and Metrology.
	This symbol means the product is certified by TUV Rheinland.
	This symbol means that product complies with Product Safety Appliance & Material Standards.
	This symbol means that the product has been approved by the Canadian Standards Association.
	This symbol means that the product complies with European Norms Electrical Certification.
	RCM Australian Compliance symbol, indicating compliance with applicable electrical safety, EMC, EME and telecommunications requirements.
	Indicates the product is a Medical Device according to EU Medical Device Regulation 2017/745.

## Safety Instructions

### Battery and Battery Charger Safety Practices

- The charger must be placed in a low traffic area to prevent it from being hit.
- When plugging / unplugging the charger's power cord, grip the plug itself. Do not pull on the cord.
- Attach any loose power cord to prevent cable wear and to avoid strangulation.
- Install the charger away from heat sources or direct sunlight.
- The charger has been made for dry surroundings. It must not be installed in bathrooms, or other humid areas.
- Only use the charger in well ventilated areas.
- The charger must not be covered or exposed to dust or water. The metal contacts must be kept clean.
- The charger is for use only with Arjo 24V lead-acid battery type, as supplied with the lift.
- The battery pack may remain connected to the charger unit when it is fully charged without damage.
- Do not short-circuit the battery pack.
- Do not crush, puncture, open, dismantle or otherwise mechanically interfere with the batteries.
- The power socket must be easily accessible. If a faulty condition occurs, switch off the power by removing the plug from the wall socket.
- Do not store battery packs at a temperature higher than 40 °C (104 °F). Should the battery pack casing crack and cause contents to come in contact with skin or clothing, rinse immediately with water. If content comes in contact with eyes, rinse immediately with plenty of water and seek medical attention.
- The "Pb" symbol appears next to the recycling and trash bin symbols on the battery pack label. This means that the battery contains lead and that it must not be disposed of, but recycled.

# Installation Instructions

## Preparation

- Mount the charger on a suitable wall close to an electrical outlet.
- Check if wall plugs or anchors are needed. No hardware is supplied with the kit.
- Use anchors, plugs and screws (4) able to bear the weight of the charger with a battery in it. Each anchor point must resist to a pull force of 110 N (25 lb). Holes in mounting plate are 5 mm (3/16 in) diameter.

## Installation

- 1) Use the charger mounting plate, as a template, to mark the four screw's location on the wall. Make sure the charger is held straight.

- 2) Drill holes (at markings) into the wall. Fit wall plugs into the wall if needed.
- 3) Hold charger to wall, fit and tighten screws, see Fig. 1.
- 4) Plug the charger using the power cord supplied with the kit.
- 5) Attach loose power cord.

**WARNING:** Precaution shall be taken, especially in home care environment, to prevent strangulation risks by routing cables properly avoiding loops.

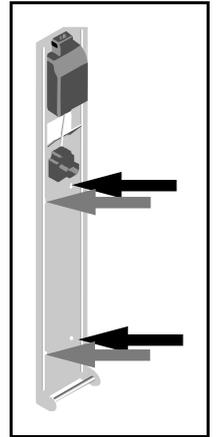


Fig. 1

## How to Use the Charger

### Battery Pack

Battery life is variable (2-3 years) and mainly depends on proper charging practices.

Recharge the batteries on a regular basis to prolong their life span. Avoid the batteries from reaching a low charge state.

Remove the battery pack from the lift when it is not used for a long period.

**NOTE:** Batteries must be charged for at least 8 hours prior to the initial use in order to reach optimal life time.

- 4) When the batteries are fully charged, remove the pack from the charger. To do so, press the red button and pull the battery pack towards you (see Fig. 3). The LED remains green to show that the charger is powered.

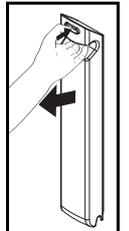


Fig. 3

**NOTE:** It is recommended to unplug the charger when it is not being used.

**CAUTION:** Refer to your local regulation for recycling and disposal of battery packs.

### Charging your Battery

- 1) Insert the battery pack into the charger unit and push it in firmly (see Fig. 2).

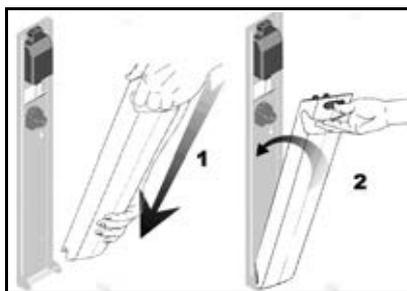


Fig. 2

**WARNING :** Hold the battery pack firmly to ensure it does not drop. Dropping can cause damages resulting in personal injuries.

- 2) Connect the charger into an AC outlet using the proper power cord.
- 3) The charger's LED turns amber during charging. The battery pack is ready for use when the LED turns green. A depleted battery pack can take up to eight hours to reach full charge state.

# Troubleshooting

The charger's LED does not light up:

- Check if the power cord is plugged into the wall outlet and into the charger.

- Check if there is power in the wall outlet.

If the LED still not lit, contact your Arjo agent to replace the charger.

# Technical Specifications

## Technical Data

PRODUCT INFORMATION	Wall Mounted Battery Charger
Charger model	NDA8200
Total weight of charger	2 kg
Compatible with battery type and part number	Arjo 24V lead-acid batteries part numbers: NDA0100-03, NDA0100-20 and NDA0200-20
Input voltage and frequency	100-240 VAC, 50-60 Hz
Power consumption	50 W (full load), ≤ 0.75 W (standby)
Output voltage and charger current	24-29 VDC, 1 A max., 24 VA
Degree of protection	IP40
Insulation class	Class II
Operating temperature / humidity range	0 °C to 40 °C (32 °F to 104 °F), 90 % RH max
Storage temperature / humidity range	-40 °C to 70 °C (-40 °F to 158 °F), 10 % RH to 95 % RH
<p>This Arjo product meets the requirements of the Low Voltage Directive 2014/35/EU, the regulations of the EMC Directive 2014/30/EU, ECO Design Directive 2009/123/EC and ROHS Directive 2011/65/EU.</p> <p>The Wall Mounted Battery Charger is compliant to IEC 60601-1 series including applicable collateral standards and national deviations. The Wall Mounted Battery Charger is also compliant to EN 60335-1:2012+A11:2014 and EN 55014-1:2006+A1:2009.</p>	
<p><b>WARNING:</b> Portable RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 30 cm to any part of the Wall Mounted Battery Charger, including cables specified by the manufacturer. Otherwise, performance degradation of this equipment could result. See "Electromagnetic Compatibility" section for more details.</p>	
<p><b>WARNING:</b> This equipment is not suitable for use in the presence of flammable anesthetic mixture with air or oxygen, or with nitrous oxide.</p>	
Safe Disposal at End of Life	
Package	Cardboard, recyclable
Electrical and electronic components	Charger systems having electrical and electronic components or an electrical cord should be disassembled and recycled per Waste of Electrical and Electronic Equipment (WEEE) or in accordance with local or national regulation.

# Technical Specifications

## Dimensions

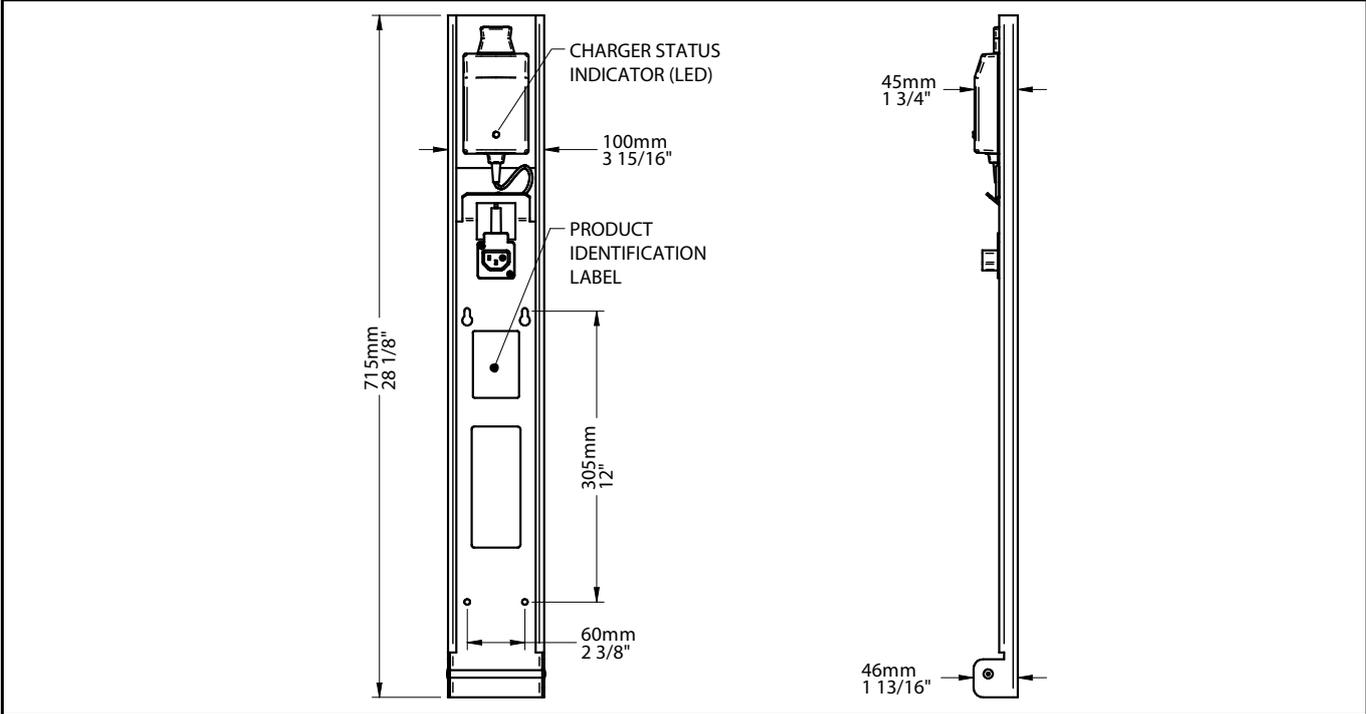


Fig. 4

# Electromagnetic Compatibility

## Electromagnetic Compliance

The Wall Mounted Battery Charger has been tested for compliance with current regulatory standards regarding its capacity to block EMI (electromagnetic interference) from external sources.

Nonetheless, some procedures can help reduce electromagnetic interferences:

- Ensure that other devices in patient-monitoring and/or life-support areas comply to accepted emissions standards.
- Maximize the distance between electro-medical devices. High-powered devices may produce EMI that can affect the charger.

For more information on how to manage the unit's RF electromagnetic environment, please consult the *AMI TIR 18-1997 - Guidance on Electromagnetic Compatibility of Medical Devices for Clinical/Biomedical Engineers*.

**WARNING : Use of this equipment adjacent to or stacked with other equipment should be avoided because it could result in improper operation. If such use is necessary, this equipment and the other equipment should be observed to verify that they are operating normally.**

**WARNING : The equipment may cause radio interference or may disrupt the operation of nearby equipment. It may be necessary to take action, such as reorienting, relocating the equipment or shielding the location.**

**WARNING : Use of accessories, cables and spare parts other than those specified or provided by Arjo could result in increased electromagnetic emissions or decreased electromagnetic immunity of this equipment and result in improper operation.**

## Electromagnetic Emissions

### Guidance and Manufacturer's Declaration - Electromagnetic Emissions - For all Equipment and Systems

The Wall Mounted Battery Charger is intended for use in the electromagnetic environment indicated below. The customer or the user of the Wall Mounted Battery Charger should assure that it is used in such an environment.

Emissions test	Compliance	Electromagnetic environment - guidance
RF emissions CISPR 11	Group 1	The Wall Mounted Battery Charger uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.  The Wall Mounted Battery Charger is suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
RF emissions CISPR 11	Class B	
Harmonic emissions IEC 61000-3-2	Class A	
Voltage fluctuations/flicker emissions IEC 61000-3-3	Complies	

# Electromagnetic Compatibility

## Electromagnetic Immunity

### Guidance and Manufacturer's Declaration - Electromagnetic Immunity - For all Equipment and Systems

The Wall Mounted Battery Charger is intended for use in electromagnetic environment specified below. The customer or the user of the Wall Mounted Battery Charger should assure that it is used in such an environment.

Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment - guidance
Electrostatic discharge (ESD) IEC 61000-4-2	±8 kV contact ±15 kV air	±8 kV contact ±15 kV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%.
Electrical fast transient/burst IEC 61000-4-4	±2 kV, AC Mains 100 kHz repetition frequency	±2 kV, AC Mains 100 kHz repetition frequency	Mains power quality should be that of a typical commercial or hospital environment.
Surge IEC 61000-4-5	±2 kV AC Mains, Lines to Ground ±1 kV, AC Mains, Line to Line	±2 kV AC Mains, Lines to Ground ±1 kV, AC Mains, Line to Line	Mains power quality should be that of a typical commercial or hospital environment.
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	0% $U_T$ ; 0,5 cycle At 0°, 45°, 90°, 135°, 180°, 225°, 270° and 315°  0% $U_T$ ; 1 cycle and 70% $U_T$ ; 25/30 cycles Single phase: at 0°  0% $U_T$ ; 250/300 cycle	0% $U_T$ ; 0,5 cycle At 0°, 45°, 90°, 135°, 180°, 225°, 270° and 315°  0% $U_T$ ; 1 cycle and 70% $U_T$ ; 25/30 cycles Single phase: at 0°  0% $U_T$ ; 250/300 cycle	Mains power quality should be that of a typical commercial or hospital environment.
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	30 A/m 50/60 Hz	30 A/m 50/60 Hz	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.

NOTE:  $U_T$  is the AC mains voltage prior to application of the test level.

# Electromagnetic Compatibility

(continued)

Guidance and Manufacturer's Declaration - Electromagnetic Immunity - For all Equipment and Systems			
Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment - guidance
Conducted RF IEC 61000-4-6	3 V outside ISM bands between 0.15-80 MHz  6 V inside ISM and amateur radio bands between 0.15-80 MHz	3 V outside ISM bands between 0.15-80 MHz  6 V inside ISM and amateur radio bands between 0.15-80 MHz	N/A
Radiated RF IEC 61000-4-3	10 V/m 80 MHz to 2.7 GHz	10 V/m 80 MHz to 2.7 GHz	N/A
Proximity fields from RF wireless communications equipment IEC 61000-4-3	380 - 390 MHz 27 V/m; PM 50%; 18 Hz  430 - 470 MHz 28 V/m; (FM $\pm 5$ kHz, 1 kHz sine) PM; 18 Hz  800 - 960 MHz 28 V/m; PM 50%; 18 Hz  1700 - 1990 MHz 28 V/m; PM 50%; 217 Hz  2400 - 2570 MHz 28 V/m; PM 50%; 217 Hz  5100 - 5800 MHz 9 V/m; PM 50%; 217 Hz  704 - 787 MHz 9 V/m; PM 50%; 217 Hz	380 - 390 MHz 27 V/m; PM 50%; 18 Hz  430 - 470 MHz 28 V/m; (FM $\pm 5$ kHz, 1 kHz sine) PM; 18 Hz  800 - 960 MHz 28 V/m; PM 50%; 18 Hz  1700 - 1990 MHz 28 V/m; PM 50%; 217 Hz  2400 - 2570 MHz 28 V/m; PM 50%; 217 Hz  5100 - 5800 MHz 9 V/m; PM 50%; 217 Hz  704 - 787 MHz 9 V/m; PM 50%; 217 Hz	N/A

AUSTRALIA  
Arjo Australia Pty Ltd  
78, Forsyth Street  
O'Connor  
AU-6163 Western Australia  
Tel: +61 89337 4111  
Free: +1 800 072 040  
Fax: + 61 89337 9077

BELGIQUE / BELGIË  
Arjo NV/SA  
Evenbroekveld 16  
BE-9420 ERPE-MERE  
Tél/Tel: +32 (0) 53 60 73 80  
Fax: +32 (0) 53 60 73 81  
E-mail: info.belgium@arjo.be

BRASIL  
Arjo Brasil Equipamentos Médicos Ltda  
Rua Marina Ciufuli Zanfelice, 329 PB02 Galpão  
- Lapa  
São Paulo – SP – Brasil  
CEP: 05040-000  
Phone: 55-11-3588-5088  
E-mail: vendas.latam@arjo.com  
E-mail: servicios.latam@arjo.com

CANADA  
Arjo Canada Inc.  
90 Matheson Boulevard West  
Suite 300  
CA-MISSISSAUGA, ON, L5R 3R3  
Tel/Tél: +1 905 238 7880  
Free: +1 800 665 4831 Institutional  
Free: +1 800 868 0441 Home Care  
Fax: +1 905 238 7881  
E-mail: info.canada@arjo.com

ČESKÁ REPUBLIKA  
Arjo Czech Republic s.r.o.  
Na Strži 1702/65  
140 00 Praha  
Czech Republic  
Phone No: +420225092307  
e-mail: info.cz@arjo.com

DANMARK  
Arjo A/S  
Vassingerødvej 52  
DK-3540 LYNGE  
Tel: +45 49 13 84 86  
Fax: +45 49 13 84 87  
E-mail:  
dk\_kundeservice@arjo.com

DEUTSCHLAND  
Arjo GmbH  
Peter-Sander-Strasse 10  
DE-55252 MAINZ-KASTEL  
Tel: +49 (0) 6134 186 0  
Fax: +49 (0) 6134 186 160  
E-mail: info-de@arjo.com

ESPAÑA  
Arjo Ibérica S.L.  
Parque Empresarial Rivas Futura, C/Marie  
Curie 5  
Edificio Alfa Planta 6 oficina 6.1-.62  
ES-28521 Rivas Vacia, MADRID  
Tel: +34 93 583 11 20  
Fax: +34 93 583 11 22  
E-mail: info.es@arjo.com

FRANCE  
Arjo SAS  
2 Avenue Alcide de Gasperi  
CS 70133  
FR-59436 RONCQ CEDEX  
Tél: +33 (0) 3 20 28 13 13  
Fax: +33 (0) 3 20 28 13 14  
E-mail: info.france@arjo.com

HONG KONG  
Arjo Hong Kong Limited  
Room 411-414, 4/F, Manhattan Centre,  
8 Kwai Cheong Road, Kwai Chung, N.T.,  
HONG KONG  
Tel: +852 2960 7600  
Fax: +852 2960 1711

ITALIA  
Arjo Italia S.p.A.  
Via Giacomo Peroni 400-402  
IT-00131 ROMA  
Tel: +39 (0) 6 87426211  
Fax: +39 (0) 6 87426222  
E-mail: Italy.promo@arjo.com

MIDDLE EAST  
Arjo Middle East FZ-LLC  
Office 908, 9th Floor,  
HQ Building, North Tower,  
Dubai Science Park,  
Al Barsha South  
P.O Box 11488, Dubai,  
United Arab Emirates  
Direct +971 487 48053  
Fax +971 487 48072  
Email: Info.ME@arjo.com

NEDERLAND  
Arjo BV  
Biezenwei 21  
4004 MB TIEL  
Postbus 6116  
4000 HC TIEL  
Tel: +31 (0) 344 64 08 00  
Fax: +31 (0) 344 64 08 85  
E-mail: info.nl@arjo.com

NEW ZEALAND  
Arjo Ltd  
34 Vestey Drive  
Mount Wellington  
NZ-AUCKLAND 1060  
Tel: +64 (0) 9 573 5344  
Free Call: 0800 000 151  
Fax: +64 (0) 9 573 5384  
E-mail: nz.info@Arjo.com

NORGE  
Arjo Norway AS  
Olaf Helsets vei 5  
N-0694 OSLO  
Tel: +47 22 08 00 50  
Faks: +47 22 08 00 51  
E-mail: no.kundeservice@arjo.com

ÖSTERREICH  
Arjo GmbH  
Lemböckgasse 49 / Stiege A / 4.OG  
A-1230 Wien  
Tel: +43 1 8 66 56  
Fax: +43 1 866 56 7000

POLSKA  
Arjo Polska Sp. z o.o.  
ul. Ks Piotra Wawrzyniaka 2  
PL-62-052 KOMORNIKI (Poznań)  
Tel: +48 61 662 15 50  
Fax: +48 61 662 15 90  
E-mail: arjo@arjo.com

PORTUGAL  
Arjo em Portugal  
MAQUET Portugal, Lda.  
(Distributor Exclusivo)  
Rua Poeta Bocage n.º 2 - 2G  
PT-1600-233 Lisboa  
Tel: +351 214 189 815  
Fax: +351 214 177 413  
E-mail: Portugal@arjo.com

SUISSE / SCHWEIZ  
Arjo AG  
Fabrikstrasse 8  
Postfach  
CH-4614 HÄGENDORF  
Tél/Tel: +41 (0) 61 337 97 77  
Fax: +41 (0) 61 311 97 42

SUOMI  
Arjo Scandinavia AB  
Riihitontuntie 7 C  
02200 Espoo  
Finland  
Puh: +358 9 6824 1260  
E-mail: Asiakaspalvelu.finland@arjo.com

SVERIGE  
Arjo International HQ  
Hans Michelsensgatan 10  
SE-211 20 MALMÖ  
Tel: +46 (0) 10 494 7760  
Fax: +46 (0) 10 494 7761  
E-mail: kundservice@arjo.com

UNITED KINGDOM  
Arjo UK and Ireland  
Houghton Hall Park  
Houghton Regis  
UK-DUNSTABLE LU5 5XF  
Tel: +44 (0) 1582 745 700  
Fax: +44 (0) 1582 745 745  
E-mail: sales.admin@arjo.com

USA  
Arjo Inc.  
2349 W Lake Street Suite 250  
US-Addison, IL 60101  
Tel: +1 630 307 2756  
Free: +1 800 323 1245 Institutional  
Free: +1 800 868 0441 Home Care  
Fax: +1 630 307 6195  
E-mail: us.info@arjo.com

JAPAN  
Arjo Japan K.K.  
東京都港区虎ノ門三丁目7番8号 ランディック  
第2虎ノ門ビル9階  
電話: +81 (0)3-6435-6401

At Arjo, we are committed to improving the everyday lives of people affected by reduced mobility and age-related health challenges. With products and solutions that ensure ergonomic patient handling, personal hygiene, disinfection, diagnostics, and the effective prevention of pressure ulcers and venous thromboembolism, we help professionals across care environments to continually raise the standard of safe and dignified care. Everything we do, we do with people in mind.



ArjoHuntleigh AB  
Hans Michelsensgatan 10  
211 20 Malmö, Sweden  
[www.arjo.com](http://www.arjo.com)

**arjo**



CE