




UNDERPADS

Dri-Flo®


Incontinence protection pads designed to provide air permeability and absorbancy which may contribute to microclimate management. Tri-dimensional matting creates air-permeable durability.

Key features


Dri-Flo Underpads are constructed using a three-layer bonding process that integrates the top, core and backsheet of the underpad. This process produces a highly absorbent, breathable underpad that allows patients to receive the full benefits of underlying microclimate management therapies.¹

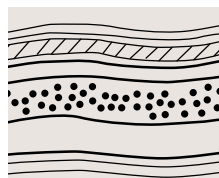

Dry

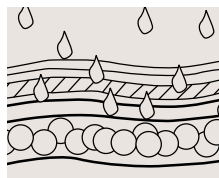

First Void

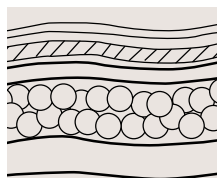

Second Void

FEATURE	BENEFIT
Soft, permeable top layer	Allows fluid to pass to the core of the underpad. Keeps moisture away from the patient's skin.
Tri-dimensional matting	Provides tensile strength and durability.
Super absorbent polymers	Draws fluid away from patient's skin and encapsulates it within the core. Absorbs secondary voids reducing associated odors.
Breathable back sheet	Allows underlying microclimate therapies to function properly. ¹

Dri-Flo Underpads:

- Provide air permeability which may aid in keeping the patient's skin dry and comfortable.
- Maintain structural integrity due to their tensile strength.
This may decrease the number of underpad changes required.

PRODUCT CODE	DESCRIPTION AND SIZE	CASE PACKAGING
7046	Large Dri-Flo Underpads 23" x 35" (58,4 cm x 88,9 cm)	6 bags of 12 pads (72 pads per case)
47001	Extra Large Dri-Flo Underpads 32" x 36" (81,3 cm x 91,4 cm)	6 bags of 8 pads (48 pads per case)

1. Arjo Data on File. Evaluation of Skin IQ: Cleveland Clinic - 100018399. September 2011.

March 2019 Only Arjo designed parts, which are designed specifically for the purpose, should be used on the equipment and products supplied by Arjo. As our policy is one of continuous development we reserve the right to modify designs and specifications without prior notice. ® and ™ are trademarks belonging to the Arjo group of companies.
© Arjo, 2019

At Arjo, we are committed to improving the everyday lives of people affected by reduced mobility and age-related health challenges. With products and solutions that ensure ergonomic patient handling, personal hygiene, disinfection, diagnostics, and the effective prevention of pressure injury and venous thromboembolism, we help professionals across care environments to continually raise the standard of safe and dignified care. Everything we do, we do with people in mind.

Arjo AB · Hans Michelsensgatan 10 · 211 20 Malmö · Sweden · +46 10 335 4500

www.arjo.com