

## HEIGHT-ADJUSTABLE BATH LIFT CHAIR

# Alenti<sup>®</sup>



Alenti is designed to enable a single caregiver to manage the full bathing cycle while maintaining a sound working posture throughout the task. Alenti fully integrates with our range of recumbent baths.

Designed with Albert, Barbara and Carl in mind



### **Facilitates efficient bathing routines**

Enables a single caregiver to collect the patient or resident from their bedside or wheelchair and transfer them directly into the tub. Alenti can also be used for showering, toileting and other hygiene tasks, reducing the number of required transfers for daily hygiene care routines.

### **Designed for easy cleaning and disinfection**

The seat is designed for easy cleaning and disinfection due to its one-piece molded construction and smooth, PUR-closed cell foam coated surfaces. The safety belt is also wipeable to reduce the risk of cross contamination.

### **Seated and supported throughout transfer**

The patient or resident can remain securely seated in Alenti as they are transferred in and out of the bath tub, thus reducing the risk of falls.

# Key features

## HEIGHT-ADJUSTABLE

Wide stroke range to suit a variety of resident heights. With built-in soft start function.

## BACK REST

Provides additional support and can be mounted on either side of the Alenti seat.

## ANATOMICALLY DESIGNED SEAT

Designed to enhance comfort and aid in correct positioning as resident seats themselves.

## SMOOTH-RUNNING CASTORS

## CENTRAL BRAKES

Activated at the press of a button.



## DUAL CONTROLS

2 sets of controls (on chassis and handset) enable caregiver to operate the device while staying close to the resident.

## FOLDABLE HAND GRIPS

Provide support throughout transfer and hygiene routine.

## RECHARGEABLE BATTERY


Two batteries and a charger are provided so you can always have a battery available while the other is charging.

## SCALE (OPTIONAL)

Highly accurate on-board electronic scale can be used as a convenient alternative for weighing residents.

## EMERGENCY RAISING & LOWERING

## SPACE REQUIREMENT

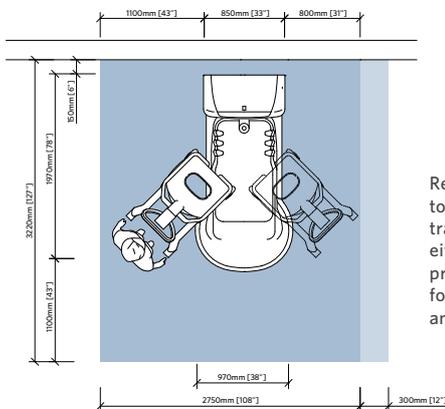

Required working area to facilitate bathing and transfer activities from either side of the tub while providing adequate access for the resident, equipment and assisting carer.

## PRODUCT INFORMATION

Safe working load	182 kg (400 lb)
Total weight (lift and scale)	72 kg (159 lb)
Total width	640 mm (25 1/4 in)
Total length	1010 mm (39 3/4 in)
Seat height, lowest position	470 mm (18 1/4 in)
Seat height, highest position	1070 (42 1/2 in)
Lifting stroke range (seat)	590 mm (23 1/4 in)
Max chassis width	640 mm (25 1/4 in)
Min chassis width between legs	478 mm (18 7/8 in)
Clearance under chassis	35 mm (1 3/8 in)
Chassis length, max (incl. brakes)	1010 mm (39 3/4 in)
Turning radius (flat stretcher position)	ø 1110 (43 3/4 in)
Supply voltage	24 V DC
Max power	200 VA
Protection class	IP X4
Protection class	IP X7 (Handset)

## STANDARD FEATURES & ACCESSORIES

- Height adjustable
- Removable back rest support
- Central brakes
- Dual controls
- Rechargeable batteries x 2
- Battery charger
- Emergency raising and lowering
- Wipeable safety belt
- Foldable hand grips

## OPTIONAL FEATURES & ACCESSORIES

- Integrated weighing scale
- Bedpan holder
- Bedpan with lid
- Leg support

For more detailed information visit [www.arjo.com](http://www.arjo.com)