

SPECIFICATION SHEET

Prioma™ Semi-Electric Couch Range



A semi-electric couch offering ergonomic and manual handling solutions

Available in 1, 2 or 3 sections with electrically operated height adjustment providing ease of use at the touch of a switch.

The 2 and 3-section models come with gas-assisted positioning of the backrest and footrest to cater to a broader range of clinical and rehabilitation needs.

Standard features

- Safe Working Load (SWL) 200 Kg
- Electrically operated variable height adjustment
- Manually operated gas assisted backrest (Prioma HC2 and Prioma HC3 only)
- Manually operated gas assisted foot rest positioning (Prioma HC3 only)
- 75mm retractable wheels
- Retractable wheel base
- Levelling foot for uneven floors
- M1 flame retardant vinyl upholstery
- Compatible with mobile hoist

Prioma HC3 couch pictured with optional breath hole.

Optional features


- ABS undercarriage cover
- Breath hole, located on the backrest (Prioma HC2 and Prioma HC3 only)

Optional accessories

- Paper roll holder
- Long side support rails - pair (optional for Prioma HC1 only)
- Short side support rails - pair (optional for Prioma HC2 and Prioma HC3 only)


Prioma HC3

3-section, adjustable height, backrest and foot positioning

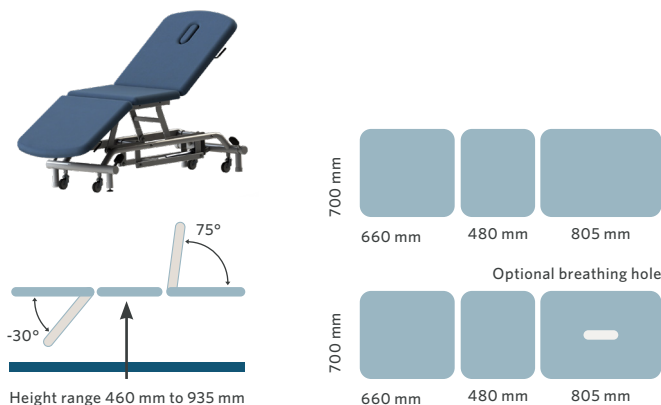


Prioma HC2

2-section, adjustable height and backrest positioning

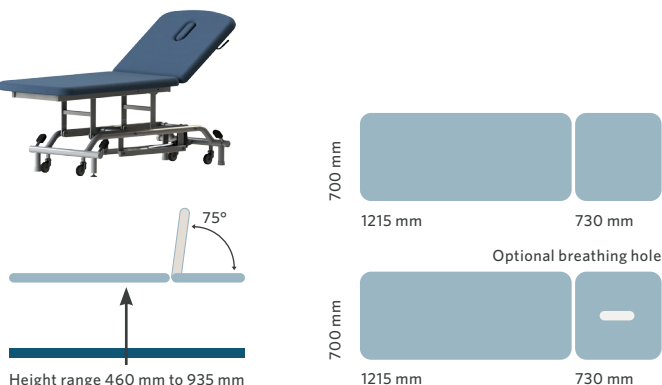

Prioma HC1

1-section, adjustable height

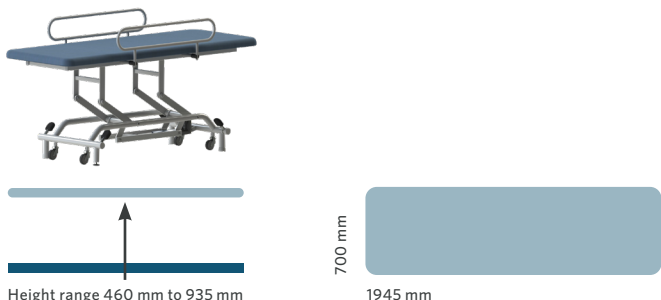

Product weight (approximate)

Safe working load (SWL)	200 Kg
Maximum patient weight (MPW)	165 Kg
Accessories weight	
Short side support rails - pair	8 Kg
Long side support rails (Prioma HC1 only) - pair	12 Kg
Paper towel rail	1.5 Kg

Electrical Data

Power in - 220V AC 50/60Hz, Max current 2.5A

Duty rating - intermittent 10% max 2 min / 18 min

Electric Shock Protection - Class 2 Type B

Liquid ingress protection - IPX4

Internally powered equipment - battery backup only - 2 x 12V series connected batteries, sealed, rechargeable lead/acid gel 1.2 amp hrs

¹ All height range measurements are from floor to top surface working height and have a tolerance of +/- 10 mm. All angular measurements have a tolerance of +/- 5°. The foot section angle depends on the height of the couch as it automatically adjusts to avoid contact with the undercarriage. The angles indicated are the maximum range and based on the couch set at a relative height. All sizes for each individual section do not include hinges and have a tolerance of +/- 10 mm.

At Arjo, we are committed to improving the everyday lives of people affected by reduced mobility and age-related health challenges. With products and solutions that ensure ergonomic patient handling, personal hygiene, disinfection, diagnostics, and the effective prevention of pressure ulcers and venous thromboembolism, we help professionals across care environments to continually raise the standard of safe and dignified care. Everything we do, we do with people in mind. www.arjo.com

Arjo AB, Hans Michelsengatan 10, SE-211 20 Malmö, Sweden. Phone: +46 (0) 10 335 45 00

www.arjo.com

arjo

Arjo.A00346.1.0.INT.EN