


**ALPHA TRANCELL  
DELUXE**  
PRESSURE REDISTRIBUTION SYSTEM

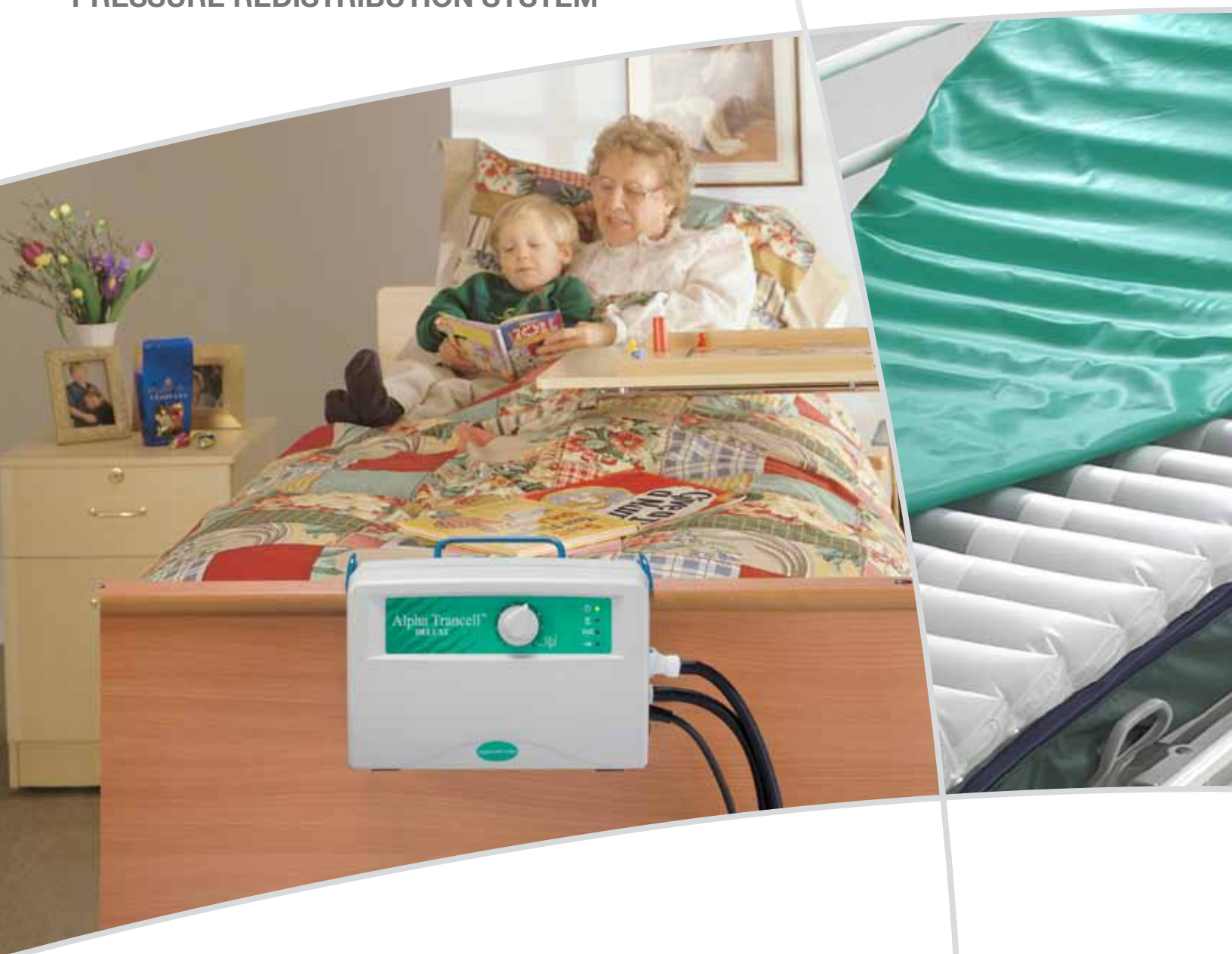


The *Alpha Trancell*® Deluxe system offers highly effective alternating pressure redistribution in a low cost, simple to use package.

Using a 10 minute cycle and a carefully calibrated pressure control on the pump, the *Alpha Trancell* Deluxe mattress replacement system is an effective tool for the prevention of pressure ulcers and in the management of superficial tissue damage.

...with people in mind

[www.ArjoHuntleigh.com](http://www.ArjoHuntleigh.com)

# KEY FEATURES

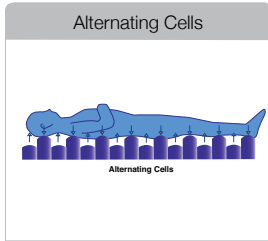

### Mattress Replacement

- 20 fully alternating modular cells
- Durable polyurethane (PU) material
- Rapid deflate function
- Transport facility
- Securing straps
- Integral foam underlay to provide full depth mattress replacement

### Pump

- Manual pressure control with easy to use numbered settings
- Lightweight and compact
- Minimal noise and vibration levels
- Tough durable casing
- On/off switch
- Power on indicator
- Inflation wait indicator
- Audible and visual power fail alarm
- Audible and visual low pressure alarm
- Integrated battery back up for alarms
- Fewer recessed areas where dirt can collect
- Hanging hooks
- Snap-fit tube connectors for secure connection
- Easy to use guide on pump

### Product information

Mattress Replacement	
Order code	ATD001DAR
Dimensions	80¼ x 34¾ x 8 1/2" (2040 x 860 x 210 mm)
Maximum Patient Weight	265 lbs (120 kg)
Cell Material	Polyurethane
Cover Material	Polyurethane coated polyester
Warranty	2 Years

Pump	
Order code	ALT03
Dimensions	9¾ x 6¼ x 4½" (248 x 160 x 116 mm)
Weight	6 lbs (2.7 kg)
Cycle	10 minutes
Power	110-120 Vac, 60 Hz, 30 VA
Electrical Safety	IEC60601-1-UK/Euro, IEC60601-1 – ROW
Warranty	2 Years



Use only ArjoHuntleigh designed parts, i.e. parts specifically designed for the purpose, on equipment and products supplied by ArjoHuntleigh. As our policy is one of continuous development, we reserve the right to modify designs and specifications without prior notice.

® and ™ are trademarks belonging to the ArjoHuntleigh group of companies. © ArjoHuntleigh, 2010